# Compact Roundness Measurement ROUNDTEST RA-10



Catalog No. E15019

Compact new roundness tester combines outstanding cost/performance ratio with full measurement capabilities



# We have a reason for proposing you introduce roundness measuring machines to your business.

Verification of geometrical tolerances, including roundness, is a must in today's quality-conscious environment.

Roundness measuring machines with the ability to perform product verification in conformity with ISO, JIS and other standards are indispensable to any quality control system that aspires to implementing high-grade quality assurance. Heightened awareness of production quality and higher quality goods will help enhance your corporate image with the buying public.

### **Roundness verification attempted** using basic measuring tools involves the following drawbacks:

- Measurement is not conducted by a radius method conforming to the standards, for which a reference axis is necessary.
- Measurement verification that meets the accuracy required by the drawings cannot be performed.
- Recorded profiles cannot be obtained.



Diameter measurement using a micrometer cannot detect an oddnumber lobing condition and resolution is marginal.




Three-point method using an indicator and V-block has better resolution but is not sensitive to common lobing

### Once roundness measuring machines are introduced into quality control:

- Reduction of nonconforming parts will translate into lower overall cost of manufacture.
- Product quality will improve and the time-to-market for new product will be reduced.
- Corporate image will be enhanced.

### **Definition of Roundness**

Roundness of a profile or contour (C) is the difference in radius (f) of two concentric circles that enclose C when the separation of these circles is a minimum, and is indicated as 'roundness xx mm' or 'roundness xx µm'.

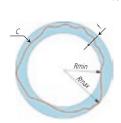





# **High-Precision Roundness Measurement**

### Simple, beginner-friendly operation

- The key layout is large and simple so is easy to view and easy to understand.
- One-shot setup recall function: Complex setups are stored in advance, ready for recall when required by one-key operation.
- Zero-setting function: The detector's level can be set to zero (0) with one single key press. This relieves the user from the chore of meticulously positioning the detector.
- The operation handles for vertical direction (Z axis) and radial direction (X axis) adjustments have been positioned on the slider for best operability.
- Because setups can only be altered in administrator mode, the machine operator can be prevented from inadvertently changing settings.

### High accuracy even though a low-end machine

Despite being a low-priced model, the turntable with air bearings offers rotational accuracy as high as (0.04+6H/10000)µm, thus assuring a precision that compares well to that of high-end models.

Large LCD panel displays measurement results and recorded profiles in an easy-to-view fashion

The built-in high-grade thermal printer prints out measurement results and recorded profiles on demand

### Compact design means small installation space

The machine calls for only a small installation space as its compact body integrates the measuring unit, electronics and printer.

### Options that further enhance usability

Use of a part setting jig exactly fitting the object being measured eliminates the need for the centering and leveling adjustments which would otherwise be required prior to measurement. An X-axis stop in the radial direction allows the detector to be positioned easily according to the object to be measured, eliminating the task of fine positioning when measurement is repeated.

\* For details on the options, see pages 3 and 8.



### Four easy steps to measurement



Clamp the workpiece to the jig.



Bring the detector into contact with the workpiece.

\*Combined use of the zero-setting function and X-axis stop (Optional) will result in securing even higher efficiency when identical workpieces are measured repetitively.



Press the [CONDITION (setup recall)] button, as needed.

\*If measurement is always conducted using the last setup, there is no need to recall this because the machine always starts up with the same settings that were effective immediately before the machine was powered down last time.



Press the [START] button.

# **Main Measuring Unit**

### **Detector**

Allows simple positioning of the workpiece due to its wide measuring range of ±1000µm.

### Part setting jig (Optional)

Can be selected to best suit the workpiece, which can be clamped/ released in a single action. High regripping accuracy eliminates the need for centering and leveling.

# High-precision air bearings

The highest accuracy in its class, (0.04+6H/10000)µm, has been achieved.

### **Built-in printer**

Prints measurement results.

### **Space-saving design**

The compact body integrating the measuring unit, electronics and printer poses no problem in installing the machine.



### **Z-axis ABS scale (Optional)**

When the ABS scale is fitted, positioning in the Z-axis (vertical) direction is performed with higher accuracy.

### X-axis stop (Optional)

Allows fast positioning of the stylus after the workpiece is clamped so that measurement can be started immediately without the need for a delicate positioning operation. This greatly increases work efficiency on batch work.

### Slider

Carries the manual operation knobs positioned together for convenient X-and Z-axis stylus position adjustment.

### Large LCD panel

Clearly displays measurement results and recorded profiles.


### Simple operation panel

Large-sized buttons allow easy recall of stored measurement setups and help prevent input errors.

## High-precision air bearings provide highly accurate measurement

Turntable axis stability is the most critical specification of a roundness measuring machine since this axis provides the datum from which the stylus deflection is measured for every type of analysis. For this reason the RA-10 is equipped with specially designed air bearings that assure high rotational accuracy to quarantee high-precision measurements.

As these bearings are inherently non-contacting they are free of any degradation arising from normal use, so the machine retains high accuracy even when used for an extended period of time.

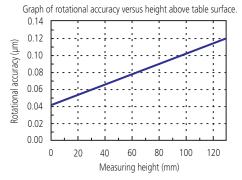



# **M**itutoyo

# Measurement results can be sent to the built-in printer or exported for external processing and storage

Measurement results and recorded profiles can be sent to the high-grade built-in thermal printer or exported via the SPC and RS-232C output functions or text file output function to USB memory.

Sample print by built-in printer

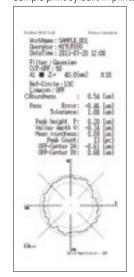

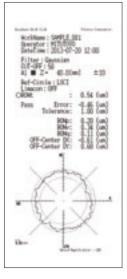

# **Control Panel**

# Measurement screen / Result screen switching

Switches between measurement screen and analysis result screen at one touch of a button.


### **Printer control**

While automatic print is available, setting can also be made to print desired results only, thus resulting in the saving of paper resources.

### **Zero Set button**

A potent tool for establishing optimum positioning of the detector.

### **Setup button**

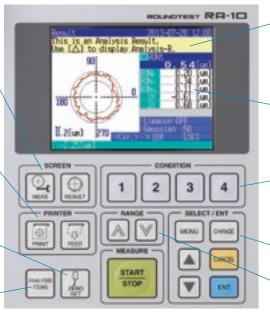

### **Supports 16 languages**

Japanese, English, German, French, Italian, Spanish, Portuguese, Korean, Traditional Chinese, Simplified Chinese, Czech, Polish, Hungarian Turkish, Swedish, Dutch

### Large LCD screen

Displays measurement results and recorded profiles in an easy-to-understand manner.

### **Setup Recall**


Frequently used measurement setups can be stored in advance, ready to be called up by one touch of a button.

### **Setup definition**

Measuring range switching

### **Simplified communication program for ROUNDTEST RA-10**

The Roundtest RA-10 has a USB interface, enabling data to be transferred to a spreadsheet or other software.

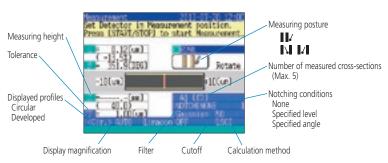

**Sample Measurement Screen** 

# Comment input Cutoff Filter Display magnification Calculation method

Sample Result Screen (Roundness)


### Useful functions help setting up prior to measurement

When a high-resolution range measurement is needed, for which accurate positioning is required, the Zero Set button allows the detector to be set at the optimum position.

The machine delivers the measurement results for a workpiece after automatically correcting for eccentricity and inclination.

### Measurement data editing function

Any part of a profile that is not to be included in the calculation can be automatically excluded from the measurement data. Therefore notches in the profile can be ignored, or data produced by scratches can be deleted while observing recorded profiles on the screen.

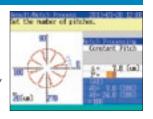

### Limaçon function compensates for eccentricity

A displacement offset between the turntable axis and that of the part under measurement results in distortion of the measured form (limaçon error) and consequentially produces an error in the calculated roundness value. The larger the eccentricity, the larger is the error in calculated roundness.

The RA-10 supports accurate measurement with a limaçon error correction function, which is provided to correct such errors arising from eccentricity.

### Notes:

- 1. The limaçon error correction is effective only when measuring a workpiece of larger diameter than that of the tip of the probe.
- If the effect obtained with the limaçon error correction function is not sufficient, use the optional alignment table (to be purchased separately) to establish precise centering of the workpiece.

# Types of Analysis with the RA-10




| Feature<br>Characteristic | Characteristic<br>Symbol    | Measurement Method | Sample Result Screen   | Explanation   |
|---------------------------|-----------------------------|--------------------|--|---|
| Form                      | Roundness                   |                    | English Analysis Brode   | Roundness (MZC) Roundness (MZC definition) of a profile is the difference in radius of two concentric circles that enclose the profile when the separation of these circles is a minimum. |
| Form                      | Flatness                    |                    | Exact In   | Flatness  Flatness of a profile is the distance between two planes enclosing the profile when this distance is a minimum.   |
| Location                  | Concentricity               |                    | the court free party benefit to complete the court of the | Concentricity  Concentricity of a profile is twice the shortest distance between the center of the profile and the datum.   |
| Location                  | Coaxiality                  |                    | the car desirate bands to the car of the car | Coaxiality  Coaxiality of the axis of a profiled surface is twice the shortest radial distance between the axis and the datum at the measured positions.                                  |
| Runout                    | Circular runout<br>(radial) |                    | Single Andrew Book and the Control of the Control o | Circular runout  Circular runout (radial) of a profile is the radial distance between two circles enclosing the profile and concentric with the datum when this distance is a minimum.    |



# **Optional Accessories**

Interchangeable Styli Unit: mm

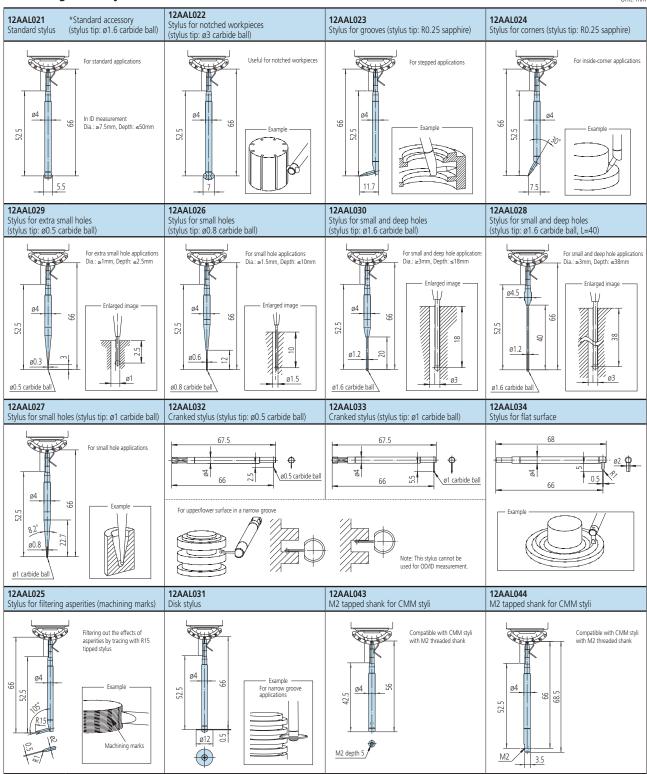

<sup>\*</sup> portion shows stylus except for the cranked stylus and stylus for flat surface.
\*( ) dimension shows a distance from the tip end of stylus or the center of tip ball to the connecting surface of detector. \*( ) dimension shows a distance from the tip end of stylus or the center of tip ball to the connecting surface of occessor.
\*Customized special interchangeable styli are available on request. Please contact any Mitutoyo office for more information.

# **Optional Accessories**

### **Various Clamping Jigs** For direct mounting on the turntable

### **Centering chuck**

When measuring a small-sized workpiece, the chuck provides good operability and the knurled ring allows the workpiece to be clamped easily.

| Order No.          | 211-052   |
|--------------------|---|
| Part holding range | O.D. (Internal jaws) 1–36mm<br>O.D. (External jaws) 25–79mm<br>I.D. (Internal jaws) 16–69mm |
| Centering error    | Within 150µm*1  |
| Mass               | 2.5kg   |

<sup>\*1:</sup> When measured with ø10mm pin gauge at measuring height of 30mm.

### V-block jig A [Semi-custom product]

The cylindrical surface of the workpiece is held against the V-block and secured with the screw-type clamp. This is a semi-custom-made product (ø10 to ø100mm) that is shipped out after adjusting the position of the V-block according to the workpiece size. This jig allows workpieces of the same size to be measured without having to center each one.

211-053: for ø50mm



### **Collet chuck**

Provides high clamping repeatability due to the use of optional precision collets. (See table below.)



<sup>1.4</sup>kg \*2: Optional collets to match the workpiece size range are required.
\*3: When measured with ø5mm pin gauge at measuring height of 30mm

### V-block jig B [Semi-custom product]

The cylindrical surface of the workpiece is held against the V-block and secured with the screw-type clamp. This is a semi-custom-made product (ø10 to ø100mm) that is shipped out after adjusting the position of the V-block according to the workpiece size. This jig allows workpieces of the same size

to be measured without having to center each one

211-054: for ø50mm

### Individual collets\*4


These collets are acquired to match the workpiece diameter range required.

| Order No.         Part Holding Range (O.D.)           12AAH402         Ø0.5-1.0mm           12AAH403         Ø1.0-1.5mm           12AAH404         Ø1.5-2.0mm           12AAH405         Ø2.0-2.5mm           12AAH406         Ø2.5-3.0mm           12AAH407         Ø3.0-3.5mm           12AAH408         Ø3.5-4.0mm           12AAH409         Ø4.0-5.0mm           12AAH410         Ø5.0-6.0mm |           |                           |
|---|-----------|---------------------------|
| 12AAH403     Ø1.0-1.5mm       12AAH404     Ø1.5-2.0mm       12AAH405     Ø2.0-2.5mm       12AAH406     Ø2.5-3.0mm       12AAH407     Ø3.0-3.5mm       12AAH408     Ø3.5-4.0mm       12AAH409     Ø4.0-5.0mm       12AAH410     Ø5.0-6.0mm   | Order No. | Part Holding Range (O.D.) |
| 12AAH404 Ø1.5–2.0mm 12AAH405 Ø2.0–2.5mm 12AAH406 Ø2.5–3.0mm 12AAH407 Ø3.0–3.5mm 12AAH408 Ø3.5–4.0mm 12AAH409 Ø4.0–5.0mm 12AAH410 Ø5.0–6.0mm   | 12AAH402  | ø0.5-1.0mm                |
| 12AAH405     Ø2.0-2.5mm       12AAH406     Ø2.5-3.0mm       12AAH407     Ø3.0-3.5mm       12AAH408     Ø3.5-4.0mm       12AAH409     Ø4.0-5.0mm       12AAH410     Ø5.0-6.0mm   | 12AAH403  | ø1.0-1.5mm                |
| 12AAH406     Ø2.5-3.0mm       12AAH407     Ø3.0-3.5mm       12AAH408     Ø3.5-4.0mm       12AAH409     Ø4.0-5.0mm       12AAH410     Ø5.0-6.0mm   | 12AAH404  | ø1.5-2.0mm                |
| 12AAH407 Ø3.0-3.5mm<br>12AAH408 Ø3.5-4.0mm<br>12AAH409 Ø4.0-5.0mm<br>12AAH410 Ø5.0-6.0mm  | 12AAH405  | ø2.0-2.5mm                |
| 12AAH408     Ø3.5-4.0mm       12AAH409     Ø4.0-5.0mm       12AAH410     Ø5.0-6.0mm   | 12AAH406  | ø2.5-3.0mm                |
| <b>12AAH409</b> Ø4.0–5.0mm<br><b>12AAH410</b> Ø5.0–6.0mm  | 12AAH407  | ø3.0-3.5mm                |
| <b>12AAH410</b> Ø5.0–6.0mm  | 12AAH408  | ø3.5-4.0mm                |
| 25.5 5.5  | 12AAH409  | ø4.0-5.0mm                |
| 42 A A LI 44 4 - C O 7 O  | 12AAH410  | ø5.0-6.0mm                |
| 12AAH411 Ø6.U-7.UMM   | 12AAH411  | ø6.0-7.0mm                |
| <b>12AAH412</b> ø7.0-8.0mm  | 12AAH412  | ø7.0-8.0mm                |
| <b>12AAH413</b> Ø8.0–9.0mm  | 12AAH413  | ø8.0-9.0mm                |
| <b>12AAH414</b> ø9.0–10.0mm   | 12AAH414  | ø9.0-10.0mm               |

<sup>\*4:</sup> A collet cannot be mounted on the turntable without a collet chuck.

### OD/ID mating jig [Semi-custom product]

These jigs are specially made to locate plain sections of a workpiece so that loading/unloading is very quick. Workpiece centering is automatically provided by just one initial centering operation on the jig, when first installed, so measurement can be started as soon as the jig is loaded with a workpiece. No clamping is used so the workpiece must be heavy enough to remain stable during measurement.

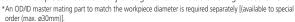



211-055: for ø10mm

### Alignment table

When installed on the turntable, this accessory enables the user to efficiently perform centering and leveling adjustments in synchronization with the adjustment Navi DAT.



| Order No.                  | 12AAH425 | 12AAH426  |
|----------------------------|----------|-----------|
| Centering adjustment range | ±3mm     | ±3mm/.12" |
| Leveling adjustment range  | ±1°      | ±1°       |
| Maximum loading            | 3kg      | 3kg       |
| Mass                       | 7kg      | 7kg       |

<sup>\*</sup>With mechanical micrometer head type (12AAH427) is available

### Options that can be installed on the alignment table

### Centering chuck (knurled ring operated)

When measuring a small-diameter workpiece, the chuck provides good operability and the knurled ring allows the workpiece to be clamped easily.



| Order No.             | 211-032  |
|-----------------------|--|
| Holding range         | O.D. with internal jaws ø1–ø36mm<br>I.D. with internal jaws ø16–ø69mm<br>O.D. with internal jaws ø25–ø79mm |
| External size (D x H) | ø118 × 41mm  |
| Mass                  | 1.2kg  |

### Microchuck

For clamping a small workpiece, 1 mm or less in diameter, that cannot be held in the centering chuck.

| Order No.             | 211-031           |
|-----------------------|-------------------|
| Holding range         | O.D.: ø0.1-ø1.5mm |
| External size (D x H) | ø107 × 48.5mm     |
| Mass                  | 0.6kg             |

### Auxiliary stage for a low-height workpiece

Used for measuring a workpiece whose diameter is 20mm or less and whose height is 20mm or less.

Order No.: 356038





### Other accessories

X-axis stop Allows the user to return the detector rapidly and easily to a fixed position on the X axis.



| Order No. | 12AAH320 |
|-----------|----------|
| Mass      | 65g      |

**SD scale for Z axis\***Scale unit for accurate positioning of the slider in the Z-axis direction (ABS scale used).



| Order No. | 12AAH318 |
|-----------|----------|
| Mass      | 450g     |

<sup>\*</sup> Shipped out attached to the RA-10 machine, or will be installed on site by Mitutoyo service personnel.

### **Vibration damping stand**



| Order No.                | 211-013                   |  |
|--------------------------|---------------------------|--|
| Vibration damping system | Diaphragm type air spring |  |
| External size            | 615 × 515 × 51mm          |  |
| Max. loading             | 150kg                     |  |


### • Simplified communication program for ROUNDTEST RA-10


The Roundtest RA-10 has a USB interface, enabling data to be transferred to a spreadsheet or other software. We also provide a program that lets you create inspection record tables using a Microsoft Excel\* macro.







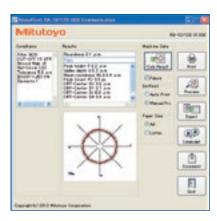

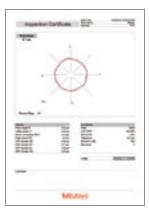

### Required environment:

- OS: Windows XP-SP3 Windows 7
- Spreadsheet software: Microsoft Excel 2010
- \*Windows OS and Microsoft Excel are products of Microsoft Corporation.

### The optional USB cable is also required.

• USB cable for **RA-10** series Order No. 12AAH490

# **Specifications**

### Main unit

| Model                                    |                                 | RA-10  |  |
|--|---------------------------------|--|--|
|  | Bearing type                    | Air bearing  |  |
|  | Rotational accuracy (radial)    | (0.04+6H/10000)µm H: Probing height (mm) JISB7451-1997   |  |
|  | Rotational accuracy (axial)     | (0.04+6X/10000)µm X: Probing radius (mm)   |  |
| Tomatable                                | Rotation speed                  | 6rpm   |  |
| Turntable                                | Effective table diameter        | ø150mm   |  |
|  | Maximum turntable loading       | 10kg   |  |
|  | Maximum probing diameter        | ø100mm   |  |
|  | Maximum workpiece diameter      | ø320mm   |  |
| Vertical column (Z axis))                | Vertical travel                 | 117mm  -Buttom position: Approx 35mm from the turntable top*2  -Top position: Approx. 152mm*1 from the turntable top*2   |  |
| 7 C. | Maximum probing height          | 152mm from the turntable top   |  |
|  | Maximum probing depth           | 100mm (minimum ID: ø30mm) using the standard stylus  |  |
| Horizontal arm (X axis)                  | Horizontal travel               | -25mm to 50mm  |  |
|  | Measuring force                 | 100mN (±30%)   |  |
| Datastas                                 | Standard stylus (12AAL021)      | Stylus tip: ø1.6mm carbide ball (Refer to page 7 for detailed information.)  |  |
| Detector                                 | Measuring range                 | ±1000µm  |  |
|  | Measuring direction             | Two directional (IN/OUT switchable)  |  |
|  | Measuring range                 | ±1000μm, ±100μm, ±10μm   |  |
|  | Magnification                   | x5, x10, x20, x50, x100, x200, x500, x1,000, x2,000, x5,000, x10,000, x20,000, x50,000, x100,000, x200,000   |  |
|  | Filter type                     | Phase corrected: 2CRPC75, 2CRPC50 Not phase corrected: 2CR75, 2CR50 Gaussion, filter OFF   |  |
|  | Cutoff value                    | 15upr, 50upr, 150upr, 500upr<br>15-150upr, 15-500upr, 50-500upr  |  |
| Electronic unit                          | Number of measuring sections    | 1-section to 5-section: Roundness, Coaxiality, Flatness 1-section to 3-section: Circular runout (radial) 2-section: Concentricity                                  |  |
|  | Reference circle for evaluation | LSCI, MZCI, MICI, MCCI   |  |
|  | Evaluation capability           | Roundness, Coaxiality, Concentricity, Flatness, Circular runout (radial)   |  |
|  | Data output                     | RS-232C VF, SPC, USB stick memory  |  |
|  | Display                         | LCD 117.2 × 88.4mm   |  |
|  | Printer                         | Thermal line printer, optional external printer  |  |
|  | Display languages               | Japanese, English, German, French, Italian, Spanish, Portuguese, Korean, Traditional Chinese, Simplified Chinese, Czech, Polish, Hungarian Turkish, Swedish, Dutch |  |
|  | Power supply                    | AC100 to 240V, 50/60Hz   |  |
|  | Power consumption               | 32–36W   |  |
| Others                                   | Air pressure                    | 0.39MPa  |  |
|  | Air consumption                 | 30L/min (minimum)  |  |
|  | Mass                            | 26kg   |  |

<sup>\*1:</sup> Top position will vary depending on any attachments installed. \*2: No attachments installed.

### **Standard accessories**


| Order No. Name                      |                                 | Quantity                           |
|-------------------------------------|---------------------------------|------------------------------------|
| 350366 Magnification adjusting film |                                 | 2 pcs.                             |
| 611755-04                           | Gauge block (35mm, JIS Grade 2) | 1 pc.                              |
| 11BAB941                            | Level                           | 1 pc.                              |
| 12AAL021                            | Standard stylus                 | 1 pc.                              |
| 12BAJ340 Printer paper*             |                                 | 2 rolls                            |
| _                                   | Receptacle                      | 1 pc.                              |
| _                                   | Hose band                       | 1 pc.                              |
| _                                   | Power cable                     | 1 pc.                              |
| _                                   | Leveling spanner                | 1 pc.                              |
| _                                   | Philips screwdriver             | 1 pc.                              |
| _                                   | Key wrench 0.9, 2 and 4         | 1 pc. (0.9), 2 pcs. (2), 1 pc. (4) |
| _                                   | Machine cover                   | 1 pc.                              |
| <ul> <li>User's manual</li> </ul>   |                                 | 1 copy                             |

<sup>\*12</sup>AAH181: Optional printer paper set (10 rolls)




# **Dimensions**


### **■** External dimensions

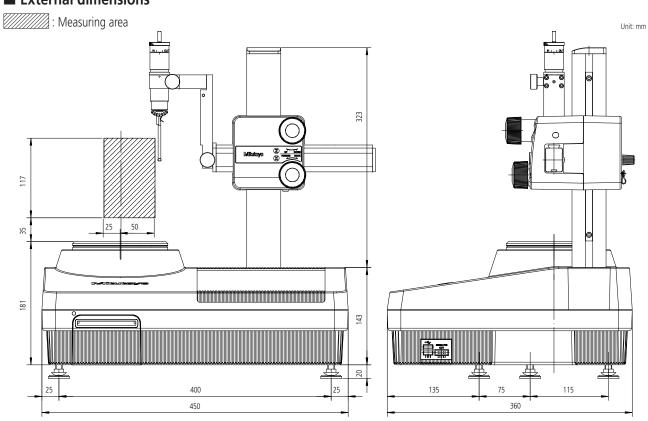

### **■** Turntable top view

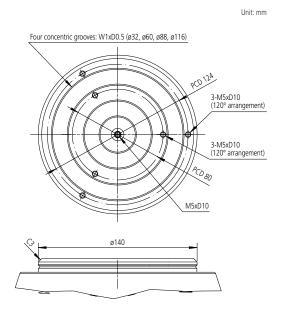

### ■ Installation floor plan

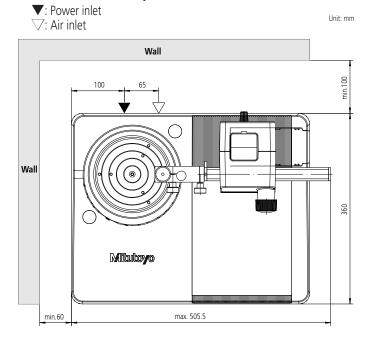



Specifications are subject to change without notice.

Export permission by the Japanese government may be required for exporting our products according to the Foreign Exchange and Foreign Trade Law. Please consult our sales office near you before you export our products or you offer technical information to a nonresident.

Coordinate Measuring Machines

Vision Measuring Systems

Form Measurement

Optical Measuring

Sensor Systems

Test Equipment and Seismometers

Digital Scale and DRO Systems

Small Tool Instruments and Data Management

### **Mitutoyo Corporation**

20-1, Sakado 1-Chome, Takatsu-ku, Kawasaki-shi, Kanagawa 213-8533, Japan T +81 (0) 44 813-8230 F +81 (0) 44 813-8231 http://www.mitutoyo.co.jp

Note: All information regarding our products, and in particular the illustrations, drawings, dimensional and performance data contained in this pamphlet, as well as other technical data are to be regarded as approximate average values. We therefore reserve the right to make changes to the corresponding designs, dimensions and weights. The stated standards, similar technical regulations, descriptions and illustrations of the products were valid at the time of printing. Only quotations submitted by ourselves may be regarded as definitive.

